

# Ambu® aScope™ 5 Broncho 4.2/2.2

## DATASHEET

aScope 5 Broncho combines the high level of manoeuvrability and imaging needed in the bronchoscopy suite with the sterility and efficiency of the single-use concept. The result? Consistently high performance, sterility, and the flexibility required for a complex and demanding workflow.

### Single-use like you've never seen it before

With the launch of the 5<sup>th</sup> generation of aScope Broncho, Ambu takes single-use bronchoscopy to even greater heights.

Building on over a decade of single-use bronchoscopy know-how, aScope 5 Broncho combines the familiar benefits of reusable scopes with the advantages of the single-use solution. It is compatible with the Ambu® aBox™ 2 and Ambu® aView™ 2 Advance display and processing units and can easily be connected to external monitors.

aScope 5 Broncho delivers consistently high performance with excellent manoeuvrability, next-level imaging and a smart, functional design. This level of performance is complemented by the sterility and efficiency of the Ambu single-use solution to give you the quality and flexibility you need for a complex and demanding workflow in the bronchoscopy suite.



## KEY BENEFITS

Completely eliminate the need for reprocessing and endoscope repairs

Start or expand bronchoscopy services with minimal capital investment

Easily navigate the bronchial segments with bending angles of 210°/210° (up/down)

Designed for smooth insertion and positioning of endotherapy instruments

Choose the right endotherapy instrument needed without fear of damaging an expensive scope

Visualize the bronchial tree clearly with the system delivering high image quality

Seamlessly integrate Ambu display and processing units to the documentation systems in use

Offer your patient a sterile scope and eliminate the risk of patient-to-patient cross-contamination

## SPECIFICATIONS

### Optical system

Field of View	120°
Direction of View	0° (forward viewing)
Depth of Field	3 - 100 mm
Illumination method	LED
Camera sensor module	CMOS

### Insertion portion

Bending section	210° up, 210° down
Insertion cord diameter	4.2 mm (0.17")
Distal tip diameter	4.4 mm (0.17")
Maximum diameter of insertion portion	Max. Ø 4.8 mm (0.19")
Minimum endotracheal tube size (inner diameter)	5.0 mm/15 Fr
Minimum DLT size	35 Fr
Working length	600 mm (23.6")
Rotary function	120° left/right
Depth marks	Every 5 cm

### Working channel

Instrument channel width	2.2 mm (0.09")
Minimum instrument channel width*	Min. Ø 2.15 mm (0.08")

### Suction connector

Connecting tube inner diameter	Ø 5.5 - 9.0 mm
--------------------------------	----------------

### Operating environment, storage, and transportation

Temperature - transportation	-10 - 55 °C (14 to 131 °F)
Temperature - operating environment	10 - 40 °C (50 to 104 °F)
Temperature - storage	10 - 25 °C (50 to 77 °F)
Relative humidity - storage	10 - 85 %
Atmospheric pressure - storage	50 - 106 kPa

### Sterilisation

Method of sterilisation	ETO
-------------------------	-----


### Other

Compatibility	Cryobiopsy probes
Latex, phthalate and PVC free	Yes

## ORDERING SPECIFICATIONS

Product code	Quantity (Ambu® aScope™ 5 Broncho 4.2/2.2)
620001000	5

\* There is no guarantee that accessories selected solely using this minimum instrument channel width will be compatible in combination.

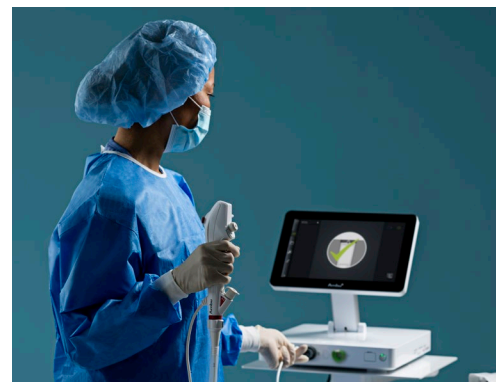

aScope 5 Broncho 4.2/2.2 works with Ambu® aBox™ 2 and aView™ 2 Advance display and processing unit and can easily be connected to external monitors.



aScope 5 Broncho 4.2/2.2 is compatible with endotherapy instruments designed to fit a working channel diameter of up to 2.2 mm.

**Ambu** FOREVER FORWARD

**Ambu Australia Pty Ltd | Ambu New Zealand Pty Ltd**

Unit 5, 2 Daydream Street, Warriewood, NSW 2102

Tel (AUS): 1300 233 118 | Tel (NZ): 0800 804 079

Email: aussales@ambu.com

www.ambuaustralia.com.au

AA23 0147 V01. Always read the Instructions for Use.

These products are only available for sale to healthcare professionals.



2797