

DATASHEET

The medium-sized Ambu BlueSensor T features conductive wet gel, offset connector and suitable adhesion, to ensure a stable signal in stress tests and short-term applications.

Key benefits

- Offset connector mitigates the impact on signal quality from external forces, like pulling of lead wires and pressure on electrodes
- Offset connector allows the user to connect and disconnect the lead wire without applying pressure to the patient or the sensor area, which reduces the risk of gel spreading and adhesive detachment from skin
- Conductive wet gel lowers skin impedance effectively and immediately
- Wet gel optimises contact between skin and electrode
- Adhesive is gentle to the skin and suitable for short-term application
- Ag/AgCl sensor ensures stable signal quality
- Visible release liner can help prevent medical accidents for both users and patients


Ambu® BlueSensor T

ECG Electrode

Single use



Wet gel

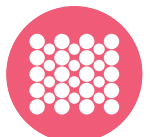

Foam backing



Silver sensor

BlueSensor

Ambu

SPECIFICATIONS

Dimensions

Electrode size (length x width in mm)	57 x 35
Skin contact size (diameter in mm)	48 x 35
Adhesive area (in mm ²)	1163

Sensor

Sensor material	Silver/silver chloride (Ag/AgCl)
Gel system	Wet gel
Gel area/measuring area (in mm ²)	254

Electrical data (ANSI/AAMI)

AC impedance - typical	450 Ω
DC offset voltage - typical	2 mV
Defibrillation overload recovery - typical	10.08 mV
Rate of change of polarisation potential - typical	-0.17 mV/s
Combined offset instability and internal noise	0.8 μV
Bias current tolerance (over 8 hours)	3.2 mV

MATERIALS

Electrode

Bio-compatible	Yes
Sponge	Polyurethane foam (PUR)
Backing material	Polyvinyl chloride foam (PVC)
Skin adhesive	Acrylate
Upper part	Polyvinyl chloride film (PVC)
A connector fitting	Polypropylene (PP)
S connector (stud)	Stainless steel
Release liner	Polyethylene film (PE)/Amorphous Polyethylene Terephthalate (APET)

AVAILABLE CONFIGURATIONS*

Reference number	Units/liner	Units/pouch	Units/inner box	Units/outer box
T-10-S/10	1	10	500	1000
T-00-A/25	1	25	500	1000
T-00-S/25	1	25	500	1000

R-00-S ■ Product ■ Wire length ■ Connector

*some configurations may not be available worldwide

CONNECTORS



S = Stud



A = 4 mm fitting

FEATURES



Printed release liner

Environment

Electrode is not made with natural rubber latex
Wire is not made with natural rubber latex
Packaging is PVC-free

X-ray & MRI

Radiolucent No
MR Conditional No

Shelf life

Opened pouches 1 month
Unopened pouches* 24 months
*from date of production

Packaging

Pouch (outer/middle/inner layers) Polyester/Aluminium/Polyethylene
Boxes Cardboard

Options

Non-wired connectors S, A

Precaution

Single use only

Ambu

Ambu Australia Pty Ltd

Unit 5, 2 Daydream Street,
Warriewood
NSW 2102
Tel: 1300 233 118
aussales@ambu.com
www.ambuaustralia.com.au

CE US: Rx only