

Ambu® aScope™ 5 Uretero

Single-use

Every case, every day

DATASHEET

aScope 5 Uretero is a single-use flexible ureteroscope solution that offers clear imaging, excellent manoeuvrability, and the flexibility needed to fit your current OR setup.

Perform procedures confidently without worry about breakage

Designed to balance robustness and manoeuvrability, the solution enables you to access difficult-to-reach areas, such as lower pole calyces. It offers excellent bending capabilities and is compatible with most common endourology instruments, such as stone removal devices, lithotripsy laser fibres, guidewires and biopsy forceps. High-resolution images make it easy to identify stones and stone fragments. **aScope 5 Uretero** is always available and never requires repairs.

A solution built for the endourology experience

The solution offers the flexibility needed to fit your existing OR setup. With a light source, a processor and a recorder all in one compact design, the compatible aBox™ 2 endoscopy system fits neatly into your existing OR setup. You can work with it independently or easily connect it to your existing visualisation tower and use it together with your external monitor. In addition, **aBox 2** integration with DICOM and PACS makes documentation smooth and efficient. You can send and retrieve imaging data directly from the endoscopy system.

The world's only ureteroscope with bioplastic material

The use of bioplastics in aScope 5 Uretero reduces the carbon footprint of the ABS plastics utilised for the production of our ureteroscope handles by 70 %. The material is sourced from by-products and waste, such as used cooking oil, which do not compete with food and agriculture production. We view the introduction of bioplastics in our endoscope handles as the first vital step towards enabling full recycling of our single-use endoscopes at scale.



KEY BENEFITS

Enables you to access difficult-to-reach areas

No cleaning or repairs needed

Flexibility to fit your existing OR setup

EHR integration makes documentation smooth and efficient

Bioplastic material lowers our carbon footprint

Ambu

SPECIFICATIONS

Optical system

Field of view	90°
Depth of field	2 - 50 mm
Illumination method	LED

Insertion portion

Bending angle*	270° up, 270° down
Insertion cord diameter	8.1 Fr / 2.7 mm (0.11")
Distal end diameter	7.9 Fr / 2.6 mm (0.10")
Maximum diameter of insertion portion	9.0 Fr / 3.0 mm (0.12")
Working length	680 mm (26.8")

Working channel

Minimum instrument channel width**	3.6 Fr / 1.2 mm (0.047") min. 1.15 mm
------------------------------------	---------------------------------------

Storage & transportation

Temperature - transportation	-10 - 55 °C
Temperature - storage***	10 - 25 °C
Humidity - transportation	10 - 95 %
Humidity - storage	10 - 85 %
Atmospheric pressure	80 - 106 kPa

Operating environment

Temperature	10 - 40 °C
Relative humidity	30 - 85 %
Atmospheric pressure	50 - 106 kPa
Altitude	≤ 2000 m

Sterilization

Method of sterilization	ETO
-------------------------	-----

Material



DEHP phthalates	Does not contain any DEHP phthalates
No allergenic substances	Not made with natural rubber latex

Dimensions

Weight	158 g
--------	-------

ORDERING INFORMATION

Product code	Product	Quantity
605001000	aScope™ 5 Uretero - Reverse Deflection	5 pcs

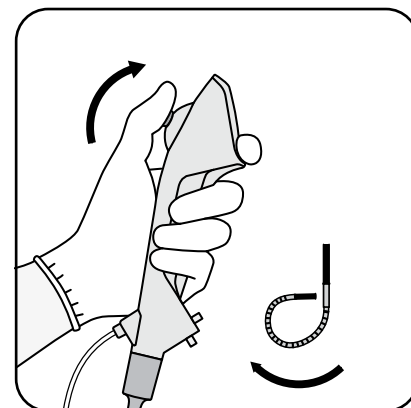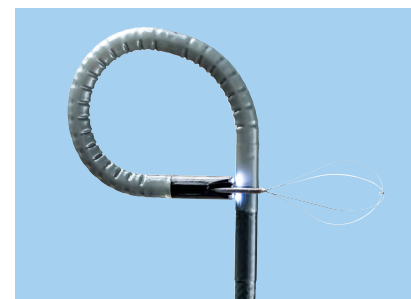

Reverse Deflection
Lever Up = Tip Down
Lever Down = Tip Up

* Please be aware that the bending angle can be affected if the insertion cord is not kept straight or endoscopic accessories have been inserted.

** There is no guarantee that endoscopic accessories selected solely using this minimum working channel width will be compatible in combination.

*** Storage under higher temperatures may impact shelf life.

Ambu FOREVER FORWARD

Ambu Australia Pty Ltd | Ambu New Zealand Pty Ltd
 Unit 5, 2 Daydream Street, Warriewood, NSW 2102
 Tel (AUS): 1300 233 118 | Tel (NZ): 0800 804 079
 Email: aussales@ambu.com Web: www.ambuaustralia.com.au

Always read the Instructions for Use.
 These products are only available for sale to healthcare professionals.

©2024 Ambu Australia Pty Ltd
 Ambu® is a registered trademark.
 aScope™ and aBox™ are trademarks of Ambu.

CE
 2797

AA24 0203 V01
 LBL-009292 - V01 - 2023/11. TCC 11567